



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
5-2-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-28-11 Cal 1 0.102 AC Cal 2 0.100 AC

Between the dates of 3-31-11 and 4-28-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


4.28-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 28, day of April, 2011.

Isabel C. Trujillo

Notary Public

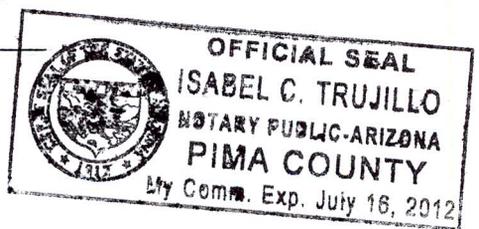

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1078  
EM  
T54  
5.2.2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 4.28-11 TIME 0837

INTOXILYZER SERIAL # 80-001199 LOCATION ODW Det HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
04/28/2011 08:38:58

2018  
CU  
BY

Standard Lot#: 719919  
Last Changed By: C. MORAND #24219

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST,SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:39:45
Diagnostic Test	Pass	08:40:12
Air Blank	0.000	08:40:39
0.100 Cal Check	0.102	08:40:58
Air Blank	0.000	08:41:27
Subject Test	0.000	08:42:17
Air Blank	ABT*	08:42:24
Air Blank	0.000	08:42:52

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
04/28/2011 08:43:13

3068  
CM  
BH

Standard Lot#: 719919  
Last Changed By: C. MORAND #24219

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:44:16
Diagnostic Test	Pass	08:44:43
Air Blank	RFI*	08:44:50
Air Blank	RFI*	08:45:05

\*RFI Detect

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
04/28/2011 08:45:39

4 of 8  
cm  
104

---

Standard Lot#: 719919  
Last Changed By: C. MORAND #24219

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:46:37
Diagnostic Test	Pass	08:47:05
Air Blank	0.000	08:47:32
0.100 Cal Check	0.100	08:47:51
Air Blank	0.000	08:48:20
Subject Test	0.111*	08:51:29
Air Blank	0.000	08:52:01
Air Blank	0.000	08:52:28
Subject Test	XXX**	08:53:34
Air Blank	0.000	08:54:07

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
04/28/2011 08:58:11

5368  
CM  
TS4

-----  
Standard Lot#: 719919  
Last Changed By: C. MORAND #24219  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:58:39
Cal Check	0.100	08:58:58
Air Blank	0.000	08:59:27

*C. Conley*

\_\_\_\_\_  
Operator's Signature

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
04/28/2011 09:00:22

Standard Lot#: 719919  
Last Changed By: C. MORAND #24219

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	09:00:57
Air Blank	0.000	09:01:25

\*Improper Sample

b of 8  
en

TS4

**Not a Successfully Completed Test Sequence**

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 8

ARIZONA DEPARTMENT OF PUBLIC SAFETY

CM

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSK  
5-2-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 4-28-11 TIME 0837  
Intoxilyzer Serial # 80-001199 Location ODW Duu Hk

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
04/28/2011 08:37:18

8569  
m

-----  
Standard Lot#: 719919  
Last Changed By: C. MORAND #24219  
-----

154

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:37:45
Cal Check	0.102	08:38:04
Air Blank	0.000	08:38:33

  
-----  
Operator's Signature