



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

By
10-29-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/11/12 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 9/13/12 and 10/11/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/18/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 18th, day of October, 2012.

Edma Escarcega

Notary Public

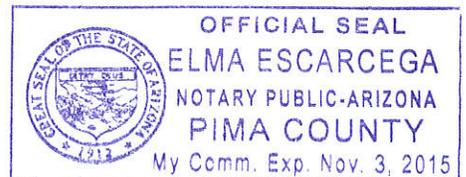

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7_{sc}

TJG
10-29-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 10/11/12 TIME 1100
INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 0.098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}
Bh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
10/11/2012 11:02:50

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:03:31
Diagnostic Test	Pass	11:04:06
Air Blank	0.000	11:04:33
0.100 Cal Check	0.099	11:04:54
Air Blank	0.000	11:05:23
Subject Test	0.000	11:05:45
Air Blank	0.000	11:06:13
Wait	ABT*	
Air Blank	0.000	11:06:48

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7 sc
BL

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
10/11/2012 11:09:19

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:10:09
Diagnostic Test	Pass	11:10:45
Air Blank	RFI*	11:10:59
Air Blank	IPS**	11:11:09

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

D4

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number: 80-001199
 Core Version: 8105.48
 10/11/2012 11:16:52

Standard Lot#: 23511100A2
 Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:17:37
Diagnostic Test	Pass	11:18:12
Air Blank	0.000	11:18:39
0.100 Cal Check	0.098	11:19:00
Air Blank	0.000	11:19:30
Subject Test	0.109*	11:22:37
Air Blank	0.000	11:23:09
Air Blank	0.000	11:23:37
Subject Test	XXX**	11:24:04
Air Blank	0.000	11:24:37

*Deficient Sample
 **Invalid Sample
 Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7_{sc}
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
10/11/2012 11:25:17

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:25:45
0.100 Cal Check	0.098	11:26:06
Air Blank	0.000	11:26:36

Preventative Maintenance
Performed -- Timer Reset

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TJ4
10.29.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/11/12 TIME 1100
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

TS4

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number: 80-001199
 Core Version: 8105.48
 10/11/2012 11:00:28

 Standard Lot#: 23511100A2
 Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:00:56
0.100 Cal Check	0.098	11:01:17
Air Blank	0.000	11:01:46

Preventative Maintenance
 Performed -- Timer Reset