



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 11.19.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/8/12 Cal 1 0.098 AC Cal 2 0.099 AC

Between the dates of 10/11/12 and 11/8/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


11/14/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14<sup>th</sup>, day of November, 2012.

Elma Escarcega

Notary Public

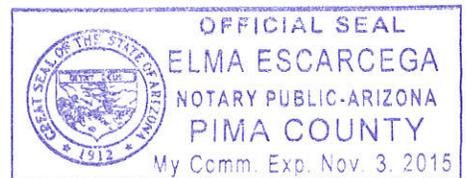

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7sc  
TS4  
11-19-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/8/12 TIME 1059

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7 sc  
BH

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
11/08/2012 11:01:43

---

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:02:26
Diagnostic Test	Pass	11:03:01
Air Blank	0.000	11:03:28
0.100 Cal Check	0.098	11:03:49
Air Blank	0.000	11:04:19
Subject Test	0.000	11:05:06
Air Blank	0.000	11:05:34
Wait	ABT*	
Air Blank	0.000	11:06:16

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7  
sc  
BH

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
11/08/2012 11:08:27

---

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	11:08:50
Air Blank	IPS**	11:08:57

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>  
134

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
11/08/2012 11:09:00

---

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:09:43
Diagnostic Test	Pass	11:10:18
Air Blank	0.000	11:10:45
0.100 Cal Check	0.099	11:11:06
Air Blank	0.000	11:11:36
Subject Test	0.273*	11:14:43
Air Blank	0.000	11:15:15
Air Blank	0.000	11:15:43
Subject Test	XXX**	11:16:06
Air Blank	0.000	11:16:39

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
11/08/2012 11:17:07

-----  
Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:17:36
0.100 Cal Check	0.099	11:17:57
Air Blank	0.000	11:18:26

Preventative Maintenance  
Performed -- Timer Reset

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TSY  
11.19.2012

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 11/8/12 TIME 1059  
Intoxilyzer Serial # 80-001199 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7<sub>sc</sub>

BH

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
11/08/2012 10:59:09

-----  
Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	α/210L	Time
Air Blank	0.000	10:59:38
0.100 Cal Check	0.098	10:59:59
Air Blank	0.000	11:00:28

Preventative Maintenance  
Performed -- Timer Reset