



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

u 8-31-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-17-2007 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 7-20-2007 and 8-17-2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.


Terry S. Gallegos

Terry S. Gallegos

8-30-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 30, day of August, 2007.

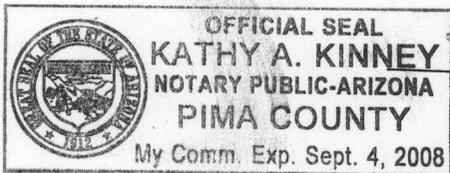

Kathy A Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

DU 8-31-07

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

1077  
TS4

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 8-17-2007 TIME 1427

INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002201  
08/17/2007

14:28:30

*Ch*  
*2 of 7*  
*TS4*

---

Standard Lot#: 706602I  
Last Changed By: J. OLIVARES #36552

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:29:26
Diagnostic Test	Pass	14:29:53
Air Blank	0.000	14:30:21
0.100 Cal Check	0.100	14:30:40
Air Blank	0.000	14:31:09
Subject Test	0.000	14:31:26
Air Blank	0.000	14:31:53
Five Minute Wait	ABT*	
Air Blank	0.000	14:32:31

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002201  
08/17/2007

OK

14:33:02

30 of 7  
124

Standard Lot#: 706602I  
Last Changed By: J. OLIVARES #36552

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:33:53
Diagnostic Test	Pass	14:34:20
Air Blank	RFI*	14:34:32
Air Blank	IPS**	14:34:46

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002201  
08/17/2007

14:37:09

CM

4077  
BY

Standard Lot#: 706602I  
Last Changed By: J. OLIVARES #36552

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:37:57
Diagnostic Test	Pass	14:38:24
Air Blank	0.000	14:38:51
0.100 Cal Check	0.099	14:39:10
Air Blank	0.000	14:39:40
Subject Test	0.135*	14:42:48
Air Blank	0.000	14:43:22
Air Blank	0.000	14:43:49
Subject Test	XXX**	14:44:06
Air Blank	0.000	14:44:40

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002201  
08/17/2007

14:45:24

*cu*

-----  
Standard Lot#: 706602I  
Last Changed By: J. OLIVARES #36552  
-----

*5077  
754*

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	14:45:51
Cal Check	0.098	14:46:10
Air Blank	0.000	14:46:40

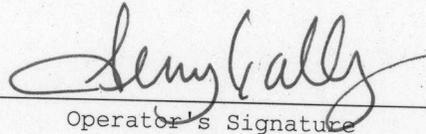
  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

CU 8-31-07  
6047  
156

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 8-17-2007 TIME 1427  
Intoxilyzer Serial # 80-002201 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002201

08/17/2007

14:27:09


-----  
Standard Lot#: 706602I

Last Changed By: J. OLIVARES #36552

-----  
QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

-----  
Test                                      g/210L                                      Time  
-----  
Air Blank                                      0.000                                      14:27:37  
Cal Check                                      0.099                                      14:27:55  
Air Blank                                      0.000                                      14:28:25

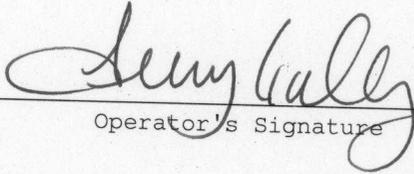  
-----  
Operator's Signature