



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
9.25.08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9.25-08 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 8.28-08 and 9.25-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


9-25-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 25th, day of September, 2008.

A.M. Hernandez

Notary Public

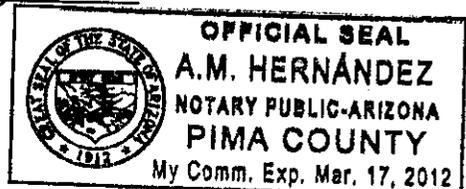

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

169
cu
Ry
9-25-08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 9-25-08 TIME 0849
INTOXILYZER SERIAL # 80-002201 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
09/25/2008

08:52:03

289 cu
PH

Standard Lot#: 735401G
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:52:54
Diagnostic Test	Pass	08:53:21
Air Blank	0.000	08:53:49
0.100 Cal Check	0.098	08:54:07
Air Blank	0.000	08:54:37
Subject Test	0.000	08:54:53
Air Blank	ABT*	08:55:00
Air Blank	0.000	08:55:28

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
09/25/2008

08:55:33

369
cu
PM

Standard Lot#: 735401G
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:56:04
Air Blank	IPS**	08:56:12

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
09/25/2008

08:57:25

Standard Lot#: 735401G
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:58:15
Diagnostic Test	Pass	08:58:42
Air Blank	0.000	08:59:09
0.100 Cal Check	0.098	08:59:28
Air Blank	0.000	08:59:58
Subject Test	0.192*	09:03:06
Air Blank	0.000	09:03:40
Air Blank	0.000	09:04:07
Subject Test	RGE**	09:04:25
Air Blank	0.000	09:05:01

*Deficient Sample

**Range Exceeded

Not a Successfully Completed Test Sequence

469
cu
104

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
09/25/2008

09:05:11

5069
cu
104

Standard Lot#: 735401G
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:05:57
Diagnostic Test	Pass	09:06:25
Air Blank	0.000	09:06:52
0.100 Cal Check	0.098	09:07:11
Air Blank	0.000	09:07:41
Subject Test	XXX*	09:08:15
Air Blank	IPS**	09:08:30

*Invalid Sample
Begin new deprivation period
**Imple

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
09/25/2008 09:09:06

bfg
on
154

Standard Lot#: 735401G
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	AMB*	09:09:34
Air Blank	0.000	09:10:03

*Ambient Fail

Clare Conley

Operator's Signature

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
09/25/2008 09:20:11

Standard Lot#: 735401G
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:20:38
Cal Check	0.097	09:20:57
Air Blank	0.000	09:21:27

7.69
cu
184

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

8049
Cu
154
9.25.08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 9-25-08 TIME 0849
Intoxilyzer Serial # 80-002201 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
09/25/2008 08:49:05

Standard Lot#: 735401G
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:49:32
Cal Check	0.097	08:49:51
Air Blank	0.000	08:50:21

Clare Conley

Operator's Signature

9569
cu
154