



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TSY  
3-17-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/12/09 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 2/12/09 and 3/12/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —

The instrument was returned to service on —

Sarah Cione

Sarah L. Cione


3/17/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17, day of March, 2009.

A.M. Hernandez

Notary Public

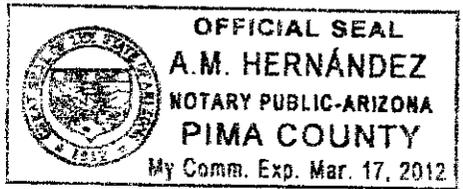

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
3/17/09

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/12/09 TIME 8:53

INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7 SC

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
03/12/2009 08:56:52

Standard Lot#: 804302I  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	α/210L	Time
Air Blank	0.000	08:57:37
Diagnostic Test	Pass	08:58:04
Air Blank	0.000	08:58:32
0.100 Cal Check	0.099	08:58:50
Air Blank	0.000	08:59:20
Subject Test	0.000	08:59:38
Air Blank	0.000	09:00:05
Five Minute Wait	ABT*	
Air Blank	0.000	09:00:39

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7sc  
TJG

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
03/12/2009 09:05:12

---

Standard Lot#: 804302I  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	09:05:56
Air Blank	IPS**	09:06:05

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
03/12/2009

09:06:21

4 of 7sc  
154

---

Standard Lot#: 804302I  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:07:34
Diagnostic Test	Pass	09:08:01
Air Blank	0.000	09:08:28
0.100 Cal Check	0.098	09:08:47
Air Blank	0.000	09:09:17
Subject Test	0.075*	09:12:25
Air Blank	0.000	09:12:58
Air Blank	0.000	09:13:25
Subject Test	XXX**	09:14:01
Air Blank	0.000	09:14:35

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/12/2009

09:16:29

-----  
Standard Lot#: 804302I  
Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:16:56
Cal Check	0.098	09:17:15
Air Blank	0.000	09:17:45

*Sarah Cione*

Operator's Signature

5 of 7 sc

B4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSY  
3-17-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/12/09 TIME 8:53  
Intoxilyzer Serial # 80-002201 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/12/2009

08:53:02

-----  
Standard Lot#: 804302I  
Last Changed By: J. OLIVARES #36552  
-----

QAS:-S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:53:30
Cal Check	0.098	08:53:48
Air Blank	0.000	08:54:18

*Sarah Cione*

Operator's Signature

7 of 7sc

TS4



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 2 sc

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TSY  
3-17-09

Date: 022109  
QAS: OLIVARES #36552  
Intoxilyzer Serial Number: 80-002201  
New Lot Number: 804302 I  
New Expiration Date: 17 FEB 2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:

Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."

- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Jose A. Olivares  
36552

2 of 2sc


154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
02/21/2009 21:11:42

-----  
Standard Lot#: 8043021  
Last Changed By: J. OLIVARES #36552

-----  
QAS: J. OLIVARES #36552  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	21:12:09
Cal Check	0.098	21:12:28
Air Blank	0.000	21:12:58

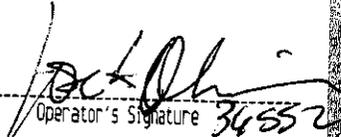-----  
  
Operator's Signature 36552