



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

104  
3-10-2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/4/10 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 2/4/10 and 3/4/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


3/4/2010

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 4<sup>th</sup>, day of March, 2010.

A.M. Hernandez

Notary Public

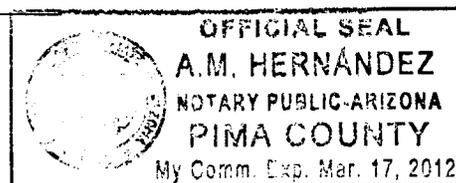

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7<sub>sc</sub>

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

BY  
3.10.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/4/10 TIME 1022

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7 sc

104

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
03/04/2010

10:23:50

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:24:33
Diagnostic Test	Pass	10:25:01
Air Blank	0.000	10:25:28
0.100 Cal Check	0.099	10:25:47
Air Blank	0.000	10:26:16
Subject Test	0.000	10:26:59
Air Blank	0.000	10:27:27
Five Minute Wait	ABT*	
Air Blank	0.000	10:28:03

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
03/04/2010

10:28:21

3 of 7 sc

TS4

---

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:29:09
Diagnostic Test	Pass	10:29:36
Air Blank	0.000	10:30:04
0.100 Cal Check	0.099	10:30:22
Air Blank	RFI*	10:30:31
Air Blank	IPS**	10:30:41

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>  
PH

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
03/04/2010

10:30:49

---

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:31:36
Diagnostic Test	Pass	10:32:04
Air Blank	0.000	10:32:31
0.100 Cal Check	0.099	10:32:50
Air Blank	0.000	10:33:20
Subject Test	0.072*	10:36:28
Air Blank	0.000	10:37:00
Air Blank	0.000	10:37:27
Subject Test	XXX**	10:37:57
Air Blank	0.000	10:38:30

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/04/2010 10:39:20

5 of 7 sc  
154

-----  
Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:39:48
Cal Check	0.098	10:40:07
Air Blank	0.000	10:40:36

  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

T54  
3.10.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/4/10 TIME 1022  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/04/2010 10:22:10

7 of 7 sc  
BU

-----  
Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:22:38
Cal Check	0.099	10:22:56
Air Blank	0.000	10:23:26

Sarah Cione  
Operator's Signature