



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
11-7-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/11 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 9/16/11 and 10/14/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 2011.

A.M. Hernandez

Notary Public

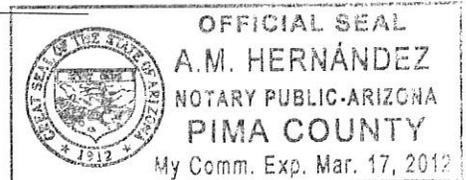

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
11-7-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/11 TIME 1413

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}

TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
10/14/2011

14:15:44

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:16:51
Diagnostic Test	Pass	14:17:26
Air Blank	0.000	14:17:53
0.100 Cal Check	0.101	14:18:14
Air Blank	0.000	14:18:44
Subject Test	0.000	14:19:12
Air Blank	0.000	14:19:39
Wait	ABT*	
Air Blank	0.000	14:20:17

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7sc

134

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
10/14/2011

14:20:26

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:21:05
Air Blank	IPS**	14:21:12

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 1027, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7_{sc}
TSH

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
10/14/2011 14:22:13

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:23:13
Diagnostic Test	Pass	14:23:48
Air Blank	0.000	14:24:15
0.100 Cal Check	0.101	14:24:36
Air Blank	0.000	14:25:06
Subject Test	0.147*	14:28:13
Air Blank	0.000	14:28:45
Air Blank	0.000	14:29:13
Subject Test	XXX**	14:29:46
Air Blank	0.000	14:30:19

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 100-102, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 7_{sc}

TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
10/14/2011 14:30:46

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:31:15
0.100 Cal Check	0.100	14:31:35
Air Blank	0.000	14:32:05

Preventative Maintenance
Performed -- Timer Reset

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
11-7-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/14/11 TIME 1413
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 7sc

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
10/14/2011 14:13:21

TS4

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:13:51
0.100 Cal Check	0.101	14:14:11
Air Blank	0.000	14:14:41

Preventative Maintenance
Performed -- Timer Reset