



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
6-27-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/21/12 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 5/24/12 and 6/21/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


6/21/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 21st, day of June, 2012.

Elma Escarcega

Notary Public

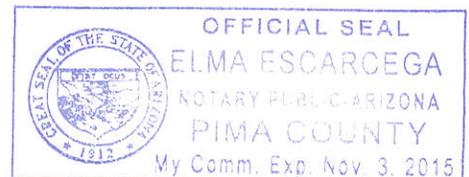

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1 of 9sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
6-27-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/21/12 TIME 1113

INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 9sc
BH

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
06/21/2012

11:19:06

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:19:53
Diagnostic Test	Pass	11:20:28
Air Blank	0.000	11:20:56
0.100 Cal Check	0.102	11:21:17
Air Blank	0.000	11:21:46
Subject Test	0.000	11:22:00
Air Blank	0.000	11:22:28
Wait	ABT*	
Air Blank	0.000	11:23:02

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 9 sc
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
06/21/2012

11:23:05

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:23:35
Air Blank	IPS**	11:23:42

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 9 sc

TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
06/21/2012 11:23:56

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:24:37
Diagnostic Test	Pass	11:25:12
Air Blank	0.000	11:25:39
0.100 Cal Check	0.101	11:26:00
Air Blank	0.000	11:26:29
Subject Test	0.117*	11:29:37
Air Blank	0.000	11:30:09
Air Blank	0.000	11:30:37
Subject Test	XXX**	11:31:15
Air Blank	0.000	11:31:48

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 9_{sc}
Bh

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
06/21/2012 11:32:11

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:32:40
0.100 Cal Check	0.100	11:33:01
Air Blank	0.000	11:33:30

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/21/12 TIME 1113
Intoxilyzer Serial # 80-002201 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

6 of 9 sc

TS4

6-27-2012

7 of 9 sc

Bh

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
06/21/2012 11:13:02

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:13:39
0.100 Cal Check	0.101	11:14:00
Air Blank	0.000	11:14:30

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

754
6-27-2012

Date: 6/21/12
 QAS: Sarah Ciont
 Intoxilyzer Serial Number: 80-002201
 New Lot Number: 235/1100A2
 New Expiration Date: 11/1/13

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press
 ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove ~~and sign printed record~~. Attach the record to the completed checklist.

sc

SIGNATURE Sarah Ciont

154

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
06/21/2012 11:35:57

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:36:26
0.100 Cal Check	0.101	11:36:47
Air Blank	0.000	11:37:17

Preventative Maintenance
Performed -- Timer Reset