



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TSy  
 10-29-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/12/12 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 9/14/12 and 10/12/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

10/18/12


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 18<sup>th</sup>, day of October, 20 12.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

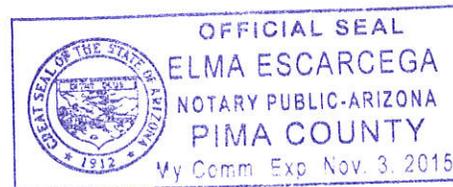

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
10-29-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/12/12 TIME 1413

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7<sub>sc</sub>  
BH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/12/2012 14:14:35

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:15:21
Diagnostic Test	Pass	14:15:56
Air Blank	0.000	14:16:23
0.100 Cal Check	0.102	14:16:44
Air Blank	0.000	14:17:13
Subject Test	0.000	14:17:30
Air Blank	0.000	14:17:58
Wait	ABT*	
Air Blank	0.000	14:18:35

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7<sub>sc</sub>  
D4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/12/2012 14:18:42

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	14:19:08
Air Blank	IPS**	14:19:16

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>  
ph

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/12/2012

14:19:20

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:20:03
Diagnostic Test	Pass	14:20:38
Air Blank	0.000	14:21:05
0.100 Cal Check	0.101	14:21:26
Air Blank	0.000	14:21:55
Subject Test	0.312*	14:25:03
Air Blank	0.000	14:25:35
Air Blank	0.000	14:26:02
Subject Test	XXX**	14:26:18
Air Blank	0.000	14:26:52

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/12/2012 14:28:17

5 of 7<sub>sc</sub>  
BH

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:28:45
0.100 Cal Check	0.101	14:29:06
Air Blank	0.000	14:29:36

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TJL  
10-29-2012

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 10/12/12 TIME 1413  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 102 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc  
15h

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/12/2012 14:13:03

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:13:33
0.100 Cal Check	0.102	14:13:53
Air Blank	0.000	14:14:23

Preventative Maintenance  
Performed -- Timer Reset