



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TSy  
1.28.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/4/13 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 12/7/12 and 1/4/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/7/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 7th, day of January, 20 13.

Edma Escarcega

Notary Public

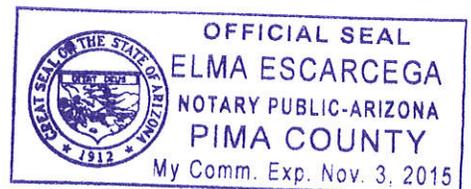

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7 SC  
TS4  
1-28-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/4/13 TIME 1336

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
01/04/2013

13:38:14

2 of 7 sc  
B4

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:39:00
Diagnostic Test	Pass	13:39:35
Air Blank	0.000	13:40:02
0.100 Cal Check	0.102	13:40:23
Air Blank	0.000	13:40:52
Subject Test	0.000	13:41:14
Air Blank	0.000	13:41:41
Wait	ABT*	
Air Blank	0.000	13:42:15

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

3 of 7 sc  
TSH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
01/04/2013 13:42:38

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:43:25
Diagnostic Test	Pass	13:44:00
Air Blank	0.000	13:44:27
0.100 Cal Check	0.102	13:44:48
Air Blank	0.000	13:45:18
Subject Test	RFI*	13:45:53
Air Blank	IPS**	13:46:01

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
01/04/2013

13:46:11

4 of 7 SC  
TS4

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:46:53
Diagnostic Test	Pass	13:47:28
Air Blank	0.000	13:47:56
0.100 Cal Check	0.102	13:48:17
Air Blank	0.000	13:48:46
Subject Test	0.047*	13:51:54
Air Blank	0.000	13:52:25
Air Blank	0.000	13:52:52
Subject Test	XXX**	13:53:14
Air Blank	0.000	13:53:48

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 SC  
TBH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
01/04/2013 13:55:24

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:55:54
0.100 Cal Check	0.101	13:56:15
Air Blank	0.000	13:56:44

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 sc  
TJG  
1-28-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 1/4/13 TIME 1336  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 102 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc  
TD4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
01/04/2013 13:36:43

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:37:12
0.100 Cal Check	0.102	13:37:33
Air Blank	0.000	13:38:02

Preventative Maintenance  
Performed -- Timer Reset