



CERTIFICATION by CUSTODIAN of RECORDS ^{W7}
Tucson Police Department Crime Laboratory ₅₋₂₅₀₁

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-26-2007 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 3-26-2007 and 4-26-2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry Gallegos


Terry S. Gallegos

5-24-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24, day of

May, 2007.

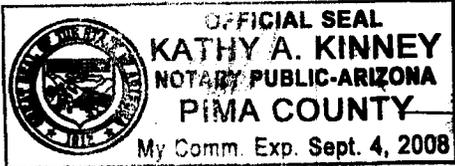

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

10 of 7
TSy 2/5
5-2007

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28047 AGENCY Tucson Police Department
DATE 4.26.2007 TIME 1513
INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
04/26/2007

15:16:13

2087
By ^{MA} 5-25-07

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:17:03
Diagnostic Test	Pass	15:17:31
Air Blank	0.000	15:17:58
0.100 Cal Check	0.100	15:18:17
Air Blank	0.000	15:18:46
Subject Test	0.000	15:19:12
Air Blank	0.000	15:19:40
Five Minute Wait	ABT*	
Air Blank	0.000	15:20:21

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
04/26/2007

15:20:37

3077
TSy
5-28-07

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:21:25
Diagnostic Test	Pass	15:21:52
Air Blank	RFI*	15:22:10
Air Blank	IPS**	15:22:22

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
04/26/2007 15:22:43

4 of 7
Bg
2ms
5/28/07

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:23:36
Diagnostic Test	Pass	15:24:03
Air Blank	0.000	15:24:31
0.100 Cal Check	0.100	15:24:50
Air Blank	0.000	15:25:19
Subject Test	0.097*	15:28:28
Air Blank	0.000	15:28:59
Air Blank	0.000	15:29:27
Subject Test	XXX**	15:29:48
Air Blank	0.000	15:30:22

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6077
TSY
5-26-07

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos AGENCY Tucson Police Department
DATE 4-26-2007 TIME 1513
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

04/26/2007

15:13:22

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

70F7
TS4
4/26/07
5:25:07

Test	g/210L	Time
Air Blank	0.000	15:13:49
Cal Check	0.099	15:14:08
Air Blank	0.000	15:14:38


Operator's Signature