



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

CU
11-23-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10-19-2007 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 9-19-2007 and 10-19-2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Terry S. Gallegos

11-21-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 21, day of November, 2007.

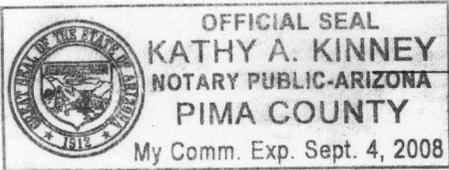

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

AM
11-23-07

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

1 of 8
Tsh

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 10-19-2007 TIME 1402

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ (code 12)

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓✓
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

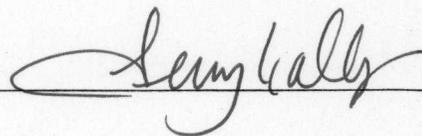SIGNATURE 

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
10/19/2007

14:04:14

cu
2078
TSy

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:05:09
Diagnostic Test	Pass	14:05:36
Air Blank	0.000	14:06:04
0.100 Cal Check	0.099	14:06:22
Air Blank	0.000	14:06:52
Subject Test	0.000	14:07:58
Air Blank	0.000	14:08:25
Five Minute Wait	ABT*	
Air Blank	0.000	14:09:00

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
10/19/2007

14:09:26

3078
754

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:10:13
Diagnostic Test	Pass	14:10:40
Air Blank	0.000	14:11:08
0.100 Cal Check	0.099	14:11:26
Air Blank	0.000	14:11:56
Subject Test	RFI*	14:12:09
Air Blank	RFI*	14:12:16

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
10/19/2007

14:12:48

4 of 8
TSh
W

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:13:49
Diagnostic Test	Pass	14:14:16
Air Blank	0.000	14:14:44
0.100 Cal Check	0.099	14:15:03
Air Blank	0.000	14:15:33
Subject Test	0.189*	14:18:41
Air Blank	0.000	14:19:13
Air Blank	0.000	14:19:41
Subject Test	XXX**	14:20:08
Air Blank	IPS***	14:20:19

*Deficient Sample

**Invalid Sample

Begin new deprivation period

***Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

10/19/2007

14:21:45

cu

Standard Lot#: 719801I

Last Changed By: T. GALLEGOS #28847

SofB
TS4


QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test g/210L Time

Air Blank AMB* 14:22:13
Air Blank 0.000 14:22:41

*Ambient Fail

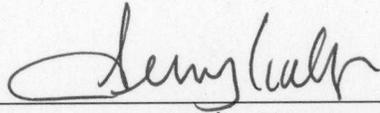

Operator's Signature

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

10/19/2007

14:23:19


Standard Lot#: 719801I

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	14:23:46
Cal Check	0.099	14:24:05
Air Blank	0.000	14:24:35

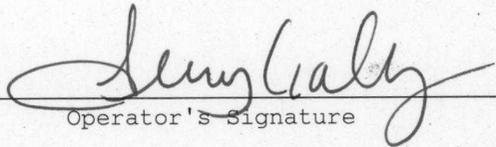
Operator's Signature

lu
bot 8
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

CM
11-23-07
7of8
TBH

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 10-19-2007 TIME 1402

Intoxilyzer Serial # 80-002261 Location DUI HQ (code 12)

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

10/19/2007

14:02:34

CM

Standard Lot#: 719801I

Last Changed By: T. GALLEGOS #28847


QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Bof 8
TSL

Test g/210L Time

Air Blank 0.000 14:03:02
Cal Check 0.100 14:03:20
Air Blank 0.000 14:03:50

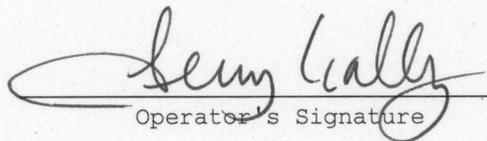

Operator's Signature