



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

By
6-10-08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/6/08 Cal 1 0.103 AC Cal 2 0.101 AC

Between the dates of 5/7/08 and 6/6/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -


Sarah Cione

Sarah L. Cione

6/9/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9, day of June, 2008.

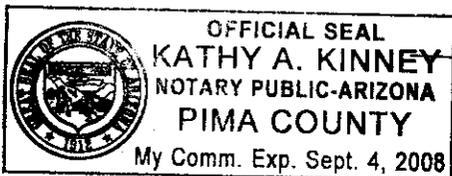

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 9 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TBY
6-10-08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/6/08 TIME 14:49
INTOXILYZER SERIAL # 80-002261 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: CMI INC
Serial Number : 80-002261
06/06/2008

14:51:10

2 of 95c

By
6.10.08

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:51:52
Diagnostic Test	Pass	14:52:20
Air Blank	0.000	14:52:47
0.100 Cal Check	0.102	14:53:06
Air Blank	0.000	14:53:36
Subject Test	0.000	14:53:55
Air Blank	0.000	14:54:23
Five Minute Wait	ABT*	
Air Blank	0.000	14:54:59

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 27-271, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: CMI INC
Serial Number : 80-002261
06/06/2008

14:55:03

3 of 95c

134

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:55:49
Diagnostic Test	Pass	14:56:17
Air Blank	RFI*	14:56:27
Air Blank	IPS**	14:56:36

*RFI Detect

**Improper Sample

4 of 95C

TS4

INTOXILYZER 8000
Location: CMI INC
Serial Number : 80-002261
06/06/2008

14:56:58

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:58:05
Diagnostic Test	Pass	14:58:33
Air Blank	0.000	14:59:00
0.100 Cal Check	0.103	14:59:19
Air Blank	0.000	14:59:48
Subject Test	0.075*	15:02:57
Air Blank	0.000	15:03:29
Air Blank	0.000	15:03:56
Subject Test	XXX**	15:04:15
Air Blank	0.000	15:04:49

*Deficient Sample
 **Invalid Sample
 Begin new deprivation period

INTOXILYZER 8000
Location: CMI INC
Serial Number: 80-002261
06/06/2008

15:05:27

5 of 9 SC
TS4

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:05:54
Cal Check	0.101	15:06:13
Air Blank	0.000	15:06:42

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 9
JC
TS4
6-10-08

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/6/08 TIME 14:49
Intoxilyzer Serial # 80-002261 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: CMI INC
Serial Number: 80-002261
06/06/2008

14:49:46

7 of 95C
TS4

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:50:13
Cal Check	0.103	14:50:32
Air Blank	0.000	14:51:01

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 SC

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 6/6/08 BY 6.10.08
 QAS: Sarah Cione
 Intoxilyzer Serial Number: 80-002261
 New Lot Number: 735401 G
 New Expiration Date: 12/28/2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 102 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: CMI INC
Serial Number: 80-002261
06/06/2008

14:48:19

9 of 95c

TS4

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:48:47
Cal Check	0.102	14:49:06
Air Blank	0.000	14:49:35

Sarah Cione

Operator's Signature