



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

SC
6/11/10

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-26-10 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 4-29-10 and 5-26-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
 The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


6-11-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 11, day of June, 2010.

A.M. Hernandez

Notary Public

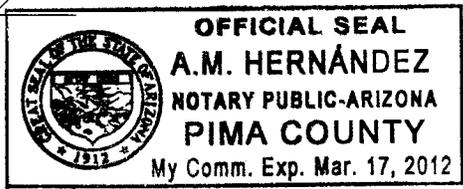

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8
EM
BY
6.1.2010

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 5-26-10 TIME 0947
INTOXILYZER SERIAL # 80-002261 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
05/26/2010

09:49:20

2578
lu
Bh

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:50:47
Diagnostic Test	Pass	09:51:15
Air Blank	0.000	09:51:42
0.100 Cal Check	0.100	09:52:01
Air Blank	0.000	09:52:30
Subject Test	0.000	09:52:59
Air Blank	ABT*	09:53:05
Air Blank	0.000	09:53:33

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
05/26/2010

09:53:39

3A8
CW
TSh

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:54:35
Diagnostic Test	Pass	09:55:03
Air Blank	RFI*	09:55:20
Air Blank	IPS**	09:55:28

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
05/26/2010

09:55:57

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:56:57
Diagnostic Test	Pass	09:57:24
Air Blank	0.000	09:57:52
0.100 Cal Check	0.100	09:58:10
Air Blank	0.000	09:58:40
Subject Test	RGE*	09:59:00
Air Blank	0.000	09:59:34

*Range Exceeded

Not a Successfully Completed Test Sequence

4588
EW
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
05/26/2010

09:59:39

5/8
lu

TS4

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:00:27
Diagnostic Test	Pass	10:00:55
Air Blank	0.000	10:01:22
0.100 Cal Check	0.100	10:01:41
Air Blank	0.000	10:02:10
Subject Test	0.093*	10:05:19
Air Blank	0.000	10:05:50
Air Blank	0.000	10:06:17
Subject Test	XXX**	10:06:47
Air Blank	0.000	10:07:20

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
05/26/2010 10:07:53

bjg
cr
TS4

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:08:20
Cal Check	0.099	10:08:39
Air Blank	0.000	10:09:09

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

758
Cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
6.1.2010

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 5-26-10 TIME 0947
Intoxilyzer Serial # 80-002261 Location ODW DUI H6

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
05/26/2010 09:47:58

848
en
B4

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:48:25
Cal Check	0.099	09:48:44
Air Blank	0.000	09:49:13

Clare Conley

Operator's Signature