



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TBY
6-25-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/24/10 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 5/26/10 and 6/24/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


6/24/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 24, day of June, 2010.

A.M. Hernandez

Notary Public

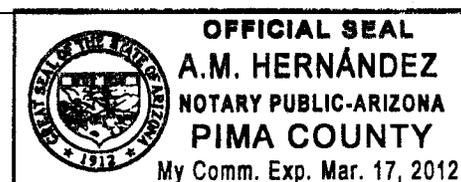

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 8 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
6-25-2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/24/10 TIME 917

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
06/24/2010

09:18:49

TS4

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:20:04
Diagnostic Test	Pass	09:20:31
Air Blank	0.000	09:20:58
0.100 Cal Check	0.100	09:21:17
Air Blank	0.000	09:21:47
Subject Test	0.000	09:22:47
Air Blank	0.000	09:23:15
Five Minute Wait	ABT*	
Air Blank	0.000	09:23:57

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8 JC

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
06/24/2010

09:26:36

104

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:27:13
Air Blank	IPS**	09:27:21

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

104

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
06/24/2010

09:27:35

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:28:17
Diagnostic Test	Fail*	09:28:21
Air Blank	0.000	09:28:49

* Analytical Stability Test

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
06/24/2010

09:29:21

TS4

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:30:27
Diagnostic Test	Pass	09:30:54
Air Blank	0.000	09:31:22
0.100 Cal Check	0.100	09:31:41
Air Blank	0.000	09:32:10
Subject Test	0.047*	09:35:18
Air Blank	0.000	09:35:49
Air Blank	0.000	09:36:17
Subject Test	XXX**	09:37:32
Air Blank	0.000	09:38:04

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

6 of 8 sc

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
06/24/2010 09:38:28

TSM

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:38:55
Cal Check	0.100	09:39:14
Air Blank	0.000	09:39:43

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

T54
6-25-2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/24/10 TIME 917
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

8 of 8 sc

D4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
06/24/2010 09:17:24

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:17:51
Cal Check	0.100	09:18:10
Air Blank	0.000	09:18:39

Sarah Cione
Operator's Signature