



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

754  
11.15.2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/11/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 10/14/10 and 11/11/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


11/11/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11 day of November, 2010.

A.M. Hernandez

Notary Public

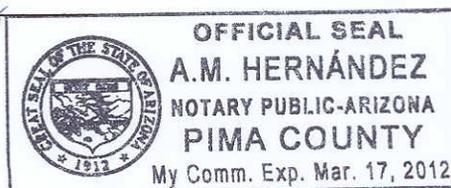

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 11 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TBY  
11.15.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 11/11/10 TIME 9:34  
INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 11 sc

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
11/11/2010 09:35:45

136

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:36:31
Diagnostic Test	Fail*	09:36:34
Air Blank	0.000	09:37:03

\* Voltage/Current Test

*- had a lot of intoxes plugged  
into outlet*

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
11/11/2010

09:39:01

3 of 11 sc

134

---

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:39:45
Diagnostic Test	Pass	09:40:12
Air Blank	0.000	09:40:40
0.100 Cal Check	0.100	09:40:58
Air Blank	0.000	09:41:28
Subject Test	0.000	09:41:43
Air Blank	ABT*	09:41:50
Air Blank	0.000	09:42:19

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

4 of 11 sc

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
11/11/2010

09:42:23

BY

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:43:09
Diagnostic Test	Pass	09:43:36
Air Blank	RFI*	09:43:44
Air Blank	IPS**	09:43:52

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

5 of 11 sc

ADH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
11/11/2010 09:44:10

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	AMB*	09:44:52
Air Blank	0.000	09:45:21

\*Ambient Fail

- had breath spray in mouth when  
improper sample was performed  
just prior to this - possible carryover

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
11/11/2010

09:46:04

6 of 11 sc  
BY

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:46:47
Diagnostic Test	Pass	09:47:14
Air Blank	0.000	09:47:42
0.100 Cal Check	0.101	09:48:00
Air Blank	0.000	09:48:30
Subject Test	0.049*	09:51:38
Air Blank	0.000	09:52:09
Air Blank	0.000	09:52:36
Subject Test	XXX**	09:53:00
Air Blank	0.000	09:53:33

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

7 of 11 sc

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
11/11/2010 09:54:08

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

TS4

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:54:36
Cal Check	0.100	09:54:55
Air Blank	0.000	09:55:24

Sarah Cione  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

8 of 11 SC

TS4  
11.15.2010

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 11/11/10 TIME 9:34  
Intoxilyzer Serial # 80-002261 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
11/11/2010 09:34:12

9 of 11 SC

SM

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:34:39
Cal Check	0.101	09:34:58
Air Blank	0.000	09:35:27

Sarah Cione  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

10 of 11 SC

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TS4  
11-15-2010

Date: 11/11/10

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002261

New Lot Number: 719919

New Expiration Date: 2/2/2012

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
  - Enter Target Value: 0.100
  - Enter Control Lot #: From the tank label
  - Enter Exp Date: From the tank label (mm/dd/yyyy)
  - Edit info Y/N? Y or N
  - Update Std Info? Y ENTER
  - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0.100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

11 of 11 SC

124

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
11/11/2010 09:56:49

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:57:17
Cal Check	0.100	09:57:35
Air Blank	0.000	09:58:05

  
-----  
Operator's Signature