



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY  
9-28-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/15/11 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 8/18/11 and 9/15/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


9/22/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of September, 2011.

A.M. Hernandez

Notary Public

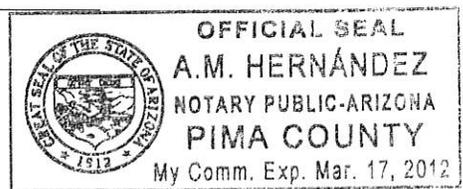

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

BH  
9.28.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/15/11 TIME 723

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Sarah Cione

2 of 7<sub>5C</sub>

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
09/15/2011 07:24:39

Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517

154

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:25:23
Diagnostic Test	Pass	07:25:50
Air Blank	0.000	07:26:17
0.100 Cal Check	0.097	07:26:36
Air Blank	0.000	07:27:06
Subject Test	0.000	07:27:21
Air Blank	0.000	07:27:48
Five Minute Wait	ABT*	
Air Blank	0.000	07:28:34

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7<sub>SC</sub>

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
09/15/2011 07:37:32

154

---

Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	07:37:58
Air Blank	IPS**	07:38:06

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>50</sub>

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
09/15/2011 07:39:33

TS4

Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:40:39
Diagnostic Test	Pass	07:41:07
Air Blank	0.000	07:41:34
0.100 Cal Check	0.097	07:41:53
Air Blank	0.000	07:42:22
Subject Test	0.103*	07:45:31
Air Blank	0.000	07:46:03
Air Blank	0.000	07:46:30
Subject Test	XXX**	07:47:00
Air Blank	0.000	07:47:32

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 7 SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
09/15/2011 07:51:38

-----  
Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517  
-----

BY

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:52:06
Cal Check	0.096	07:52:24
Air Blank	0.000	07:52:54

Sarah Cione  
Operator's Signature

6 of 7 sc

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSY  
9-28-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 9/15/11 TIME 723  
Intoxilyzer Serial # 80-002261 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7  
SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
09/15/2011 07:23:07

BY

-----  
Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:23:35
Cal Check	0.097	07:23:53
Air Blank	0.000	07:24:23



Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 2 sc

INTOXILYZER SERVICE RECORD

TS4  
9.28.2011

Date: 9/6/11

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002261

Reason For Service: officer Henemich notified crime lab that intox was not going into "Ready Mode"

Date Out of Service: N/A

Calibration Check completed prior to removal from service

Date: N/A

SQAP completed prior to removal from service

Date: N/A

Corrective Action:  Sent to CMI, Inc.

Looked at diagnostics to determine problem. Breath hose not heating up. Replaced breath hose.

159

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
09/06/2011 15:35:02

-----  
Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517  
-----

OAS: S. CIONE #50616  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:35:29
Cal Check	0.096	15:35:48
Air Blank	0.000	15:36:17

-----  
*Sarah Cione*  
Operator's Signature  
-----