



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 1-7-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/6/12 Cal 1 0.101 AC Cal 2 0.102 AC

Between the dates of 11/9/12 and 12/6/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


12/17/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17<sup>th</sup>, day of December, 2012.

Elma Escarcega

Notary Public

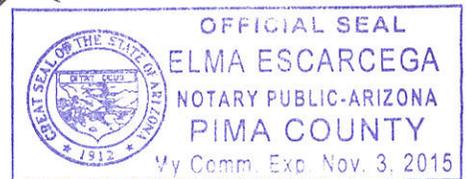

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7<sub>sc</sub>  
TS4  
1.7.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 12/6/12 TIME 1036

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

( ✓ ) 1. Clock time check.

( ✓ ) 2. Date check.

OPERATIONAL TESTS

( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.

( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 102 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7 sc  
p4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
12/06/2012

10:38:23

---

Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:39:07
Diagnostic Test	Pass	10:39:42
Air Blank	0.000	10:40:10
0.100 Cal Check	0.102	10:40:31
Air Blank	0.000	10:41:00
Subject Test	0.000	10:41:15
Air Blank	0.000	10:41:42
Wait	ABT*	
Air Blank	0.000	10:42:42

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

3 of 7<sub>sc</sub>  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
12/06/2012 10:53:32

---

Standard Lot#: 13/12100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:54:14
Diagnostic Test	Pass	10:54:49
Air Blank	0.000	10:55:17
0.100 Cal Check	0.101	10:55:38
Air Blank	0.000	10:56:07
Subject Test	RFI*	10:57:21
Air Blank	IPS**	10:57:29

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7<sub>sc</sub>

B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
12/06/2012 10:42:55

Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:43:41
Diagnostic Test	Pass	10:44:16
Air Blank	0.000	10:44:43
0.100 Cal Check	0.102	10:45:04
Air Blank	0.000	10:45:33
Subject Test	0.100*	10:48:41
Air Blank	0.000	10:49:12
Air Blank	0.000	10:49:40
Subject Test	XXX**	10:50:07
Air Blank	0.000	10:50:40

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

### Not a Successfully Completed Test Sequence

5 of 7<sub>sc</sub>  
D4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
12/06/2012 11:01:27

-----  
Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:01:55
0.100 Cal Check	0.102	11:02:16
Air Blank	0.000	11:02:45

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

6 of 7 sc  
TD4  
1-7-2013

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 12/6/12 TIME 1036  
Intoxilyzer Serial # 80-002261 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 SC  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
12/06/2012 10:36:42

-----  
Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:37:11
0.100 Cal Check	0.101	10:37:32
Air Blank	0.000	10:38:01

Preventative Maintenance  
Performed -- Timer Reset