



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

By
 1-28-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/3/13 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 12/6/12 and 1/3/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/7/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 7th, day of January, 2013.

Elma Escarcega

Notary Public

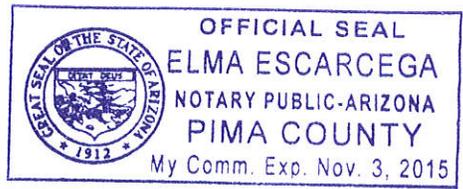

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8 sc
TS4
1-28-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/3/13 TIME 1048

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 8 sc
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/03/2013 10:51:13

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:52:09 |
| Diagnostic Test | Pass | 10:52:44 |
| Air Blank | 0.000 | 10:53:11 |
| 0.100 Cal Check | 0.102 | 10:53:32 |
| Air Blank | 0.000 | 10:54:01 |
| Subject Test | 0.000 | 10:54:41 |
| Air Blank | 0.000 | 10:55:08 |
| Wait | ABT* | |
| Air Blank | 0.000 | 10:55:48 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8 sc
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/03/2013 10:55:57

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:57:02 |
| Diagnostic Test | Pass | 10:57:37 |
| Air Blank | 0.000 | 10:58:05 |
| 0.100 Cal Check | 0.102 | 10:58:25 |
| Air Blank | 0.000 | 10:58:55 |
| Subject Test | RFI* | 10:59:09 |
| Air Blank | IPS** | 10:59:16 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 8 sc
TBH

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/03/2013 11:04:07

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 11:04:35 |
| Air Blank | 0.000 | 11:05:03 |

*RFI Detect

- opened a text message on cell phone
next to intox

Not a Successfully Completed Test Sequence

5 of 8 SC
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/03/2013

11:05:06

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F weight: iii

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:05:50 |
| Diagnostic Test | Pass | 11:06:25 |
| Air Blank | 0.000 | 11:06:53 |
| 0.100 Cal Check | 0.102 | 11:07:14 |
| Air Blank | 0.000 | 11:07:43 |
| Subject Test | 0.091* | 11:10:50 |
| Air Blank | 0.000 | 11:11:22 |
| Air Blank | 0.000 | 11:11:49 |
| Subject Test | XXX** | 11:12:10 |
| Air Blank | 0.000 | 11:12:42 |

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/03/2013 11:13:11

6 of 8 sc
TS4

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:13:39 |
| 0.100 Cal Check | 0.101 | 11:14:00 |
| Air Blank | 0.000 | 11:14:30 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8 sc
TB4
1-28-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/3/13 TIME 1048
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/03/2013 10:48:29

8 of 8 sc
TS4

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:48:57 |
| 0.100 Cal Check | 0.102 | 10:49:18 |
| Air Blank | 0.000 | 10:49:48 |

Preventative Maintenance
Performed -- Timer Reset