



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

me

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-26-2007 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 1-26-2007 and 2-26-2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos
Terry S. Gallegos

March 12, 2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 12, day of March, 2007.

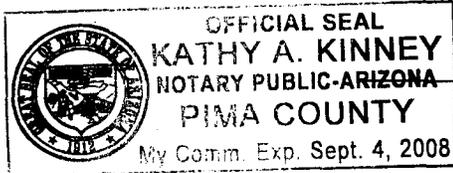

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
TSY

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 2-26-2007 TIME 1651
INTOXILYZER SERIAL # 80-002262 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002262
02/26/2007 16:51:41

2 of 7
TSY mt

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:52:33
Diagnostic Test	Pass	16:53:01
Air Blank	0.000	16:53:28
0.100 Cal Check	0.100	16:53:47
Air Blank	0.000	16:54:16
Subject Test	0.000	16:55:07
Air Blank	0.000	16:55:35
Five Minute Wait	ABT*	
Air Blank	0.000	16:56:19

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002262
02/26/2007 16:56:25

3077
By MB

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:57:20
Diagnostic Test	Pass	16:57:47
Air Blank	RFI*	16:58:04
Air Blank	IPS**	16:58:14

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002262
02/26/2007

16:58:19

4077
PH MP

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:59:11
Diagnostic Test	Pass	16:59:38
Air Blank	0.000	17:00:05
0.100 Cal Check	0.100	17:00:24
Air Blank	0.000	17:00:54
Subject Test	0.194*	17:04:02
Air Blank	0.000	17:04:34
Air Blank	0.000	17:05:01
Subject Test	XXX**	17:05:27
Air Blank	0.000	17:06:01

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002262

02/26/2007

17:06:28

5 of 7
Tg 7/17

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:06:56
Cal Check	0.100	17:07:15
Air Blank	0.000	17:07:44



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077
TSyTM

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 2-26-2007 TIME 1651
Intoxilyzer Serial # 80-00 2262 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002262

02/26/2007

16:50:15

70F7
MB
TS4

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:50:43
Cal Check	0.100	16:51:01
Air Blank	0.000	16:51:31


Operator's Signature