



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

DU
9.2.09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8.20.2009 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 7.23.2009 and 8.20.2009 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


9.2.2009

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 2nd, day of September, 20 09.

A.M. Hernández

Notary Public

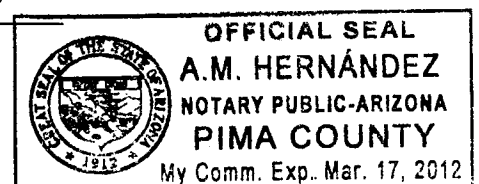

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

On 9-2-09
1 of 7
TS4

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 8.20-2009 TIME 1103

INTOXILYZER SERIAL # 80-002262 LOCATION DUT HQ

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓✓
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓
Deficient sample printed. ✓
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/20/2009

11:07:08

Standard Lot#: 826001I
Last Changed By: C. MORAND #24219

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:07:59
Diagnostic Test	Pass	11:08:26
Air Blank	0.000	11:08:53
0.100 Cal Check	0.099	11:09:12
Air Blank	0.000	11:09:42
Subject Test	0.000	11:10:07
Air Blank	ABT*	11:10:16
Air Blank	0.000	11:10:45

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/20/2009

11:11:48

Standard Lot#: 826001I
Last Changed By: C. MORAND #24219

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:12:59
Diagnostic Test	Pass	11:13:26
Air Blank	0.000	11:13:53
0.100 Cal Check	0.100	11:14:12
Air Blank	0.000	11:14:42
Subject Test	RFI*	11:15:08
Air Blank	IPS**	11:15:20

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/20/2009

11:15:43

Standard Lot#: 826001I
Last Changed By: C. MORAND #24219

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:16:46
Diagnostic Test	Pass	11:17:13
Air Blank	0.000	11:17:41
0.100 Cal Check	0.100	11:18:00
Air Blank	0.000	11:18:29
Subject Test	0.210*	11:21:37
Air Blank	0.000	11:22:10
Air Blank	0.000	11:22:37
Subject Test	XXX**	11:22:52
Air Blank	0.000	11:23:26

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
08/20/2009

11:26:10

Standard Lot#: 826001I
Last Changed By: C. MORAND #24219

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	11:26:37
Cal Check	0.099	11:26:56
Air Blank	0.000	11:27:25

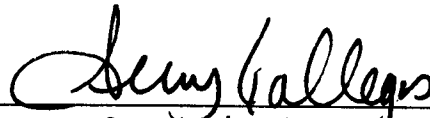

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

CH
9.2.09

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6067
TS4

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos AGENCY Tucson Police Department
DATE 8-20-2009 TIME 1103
Intoxilyzer Serial # 80-002262 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos


DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
08/20/2009 11:03:44

Standard Lot#: 826001I
Last Changed By: C. MORAND #24219

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	11:04:11
Cal Check	0.100	11:04:30
Air Blank	0.000	11:04:59

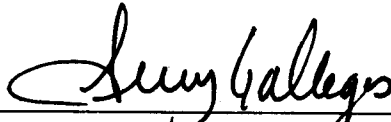

Operator's Signature