



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
8-18-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-18-11 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 7-21-11 and 8-18-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley  
Clare L. Conley


8-18-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 18<sup>th</sup>, day of August, 2011.

Susan A. Hitchye  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

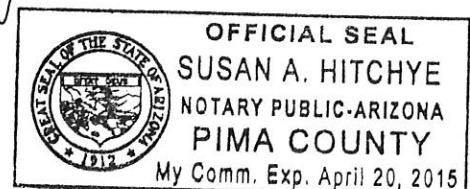

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 8-18-11 TIME 0712  
INTOXILYZER SERIAL # SD-002262 LOCATION ODW DM HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 0.00 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 1.00 AC  
Result: 0. 0.98 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
08/18/2011

07:14:11

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:15:29
Diagnostic Test	Pass	07:15:56
Air Blank	0.000	07:16:24
0.100 Cal Check	0.099	07:16:42
Air Blank	0.000	07:17:12
Subject Test	0.000	07:17:26
Air Blank	ABT*	07:17:33
Air Blank	0.000	07:18:01

\*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
08/18/2011

07:18:14

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TST, SQAP

DOB: 11/11/1111

Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:18:41
Air Blank	IPS**	07:18:48

\*RFI Detect

\*\*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
08/18/2011

07:19:11

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:19:55
Diagnostic Test	Pass	07:20:23
Air Blank	0.000	07:20:50
0.100 Cal Check	0.099	07:21:09
Air Blank	0.000	07:21:38
Subject Test	0.172*	07:24:46
Air Blank	0.000	07:25:18
Air Blank	0.000	07:25:46
Subject Test	XXX**	07:26:18
Air Blank	0.000	07:26:52

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
08/18/2011 07:27:17

-----  
Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:27:44
Cal Check	0.098	07:28:03
Air Blank	0.000	07:28:32

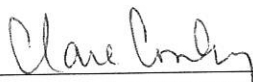
  
\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6-17  
CW

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

7/31  
8.18.2011

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 8-18-11 TIME 0712  
Intoxilyzer Serial # 80-8D2262 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
08/18/2011

07:12:48

-----  
Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:13:15
Cal Check	0.098	07:13:34
Air Blank	0.000	07:14:03

Clare Conley  
Operator's Signature