



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
2.23.2012

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/3/12 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 1/5/12 and 2/3/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/10/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 10th, day of February, 2012.

Elma Escarcega

Notary Public

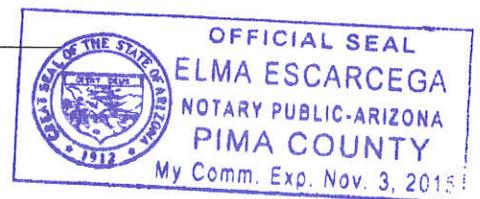

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 sc

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
2.23.2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/3/12 TIME 1431

INTOXILYZER SERIAL # 80-002262 LOCATION Vehicle #944

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

Tsh

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
02/03/2012 14:33:17

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP,SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:34:13
Diagnostic Test	Pass	14:34:42
Air Blank	0.000	14:35:09
0.100 Cal Check	0.100	14:35:28
Air Blank	0.000	14:35:58
Subject Test	0.000	14:36:14
Air Blank	0.000	14:36:41
Five Minute Wait	ABT*	
Air Blank	0.000	14:37:26

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7 sc

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
02/03/2012 14:37:34

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:38:11
Air Blank	IPS**	14:38:19

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7sc

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
02/03/2012 14:38:25

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

DAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:39:10
Diagnostic Test	Pass	14:39:39
Air Blank	0.000	14:40:06
0.100 Cal Check	0.100	14:40:25
Air Blank	0.000	14:40:54
Subject Test	0.121*	14:44:03
Air Blank	0.000	14:44:35
Air Blank	0.000	14:45:02
Subject Test	XXX**	14:45:33
Air Blank	0.000	14:46:06

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7 sc

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
02/03/2012 14:46:30

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:46:57
Cal Check	0.099	14:47:16
Air Blank	0.000	14:47:45

Sarah Cione
Operator's Signature

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

DL
2.23.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 2/3/12 TIME 1431
Intoxilyzer Serial # 80-002262 Location Vehicle # 944

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc

134

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
02/03/2012 14:31:11

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:31:38
Cal Check	0.099	14:31:57
Air Blank	0.000	14:32:27

Sarah Cione
Operator's Signature
