

TUCSON POLICE DEPARTMENT

CRIME LABORATORY *PM*

The undersigned Quality Assurance Specialist, having been qualified by the Arizona Department of Public Safety for Intoxilyzer 5000 inspection, and having performed such inspections on Intoxilyzer 5000


, Serial # 46-005003, certifies by her signature:

on 7/14/04, with a result of 0.086 AC. 0.087AC.

and on 8/11/04, with a result of 0.087 AC. 0.088AC.

this instrument was tested with a known alcohol reference standard, and always produced results within the acceptable accuracy limits (from 0.070 AC. to 0.090 AC.), as defined by R9-14-404.A.3.


No problems were reported concerning this instrument between these two dates.

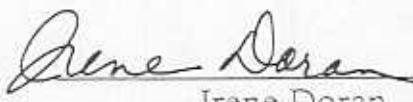


Irene Doran

ADPS QAS Permit #QAS 29586 (issue date 01/08/2004)

SUBSCRIBED AND SWORN to before me

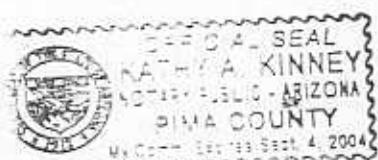

this 16, day of August, 2004.

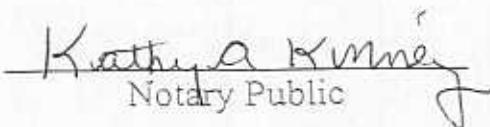

Kathy A. Kinney
Notary Public

P. log 6
JL

EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

AGENCY Tucson Police QA SPECIALIST Jennifer Combs 46160
(Print Name)
INTOXILYZER SERIAL #66-005003 LOCATION DUI HQ
DATE July 14 29 2004 TIME 1021

- Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .419.
- Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- Alcohol free subject Test result 0.000 AC
- Error recognition logic system functioning
Invalid test printed
- Proper sample recognition system
Invalid sample printed
Deficient sample printed
- Completeness of sample purge with collection tube**
- Calibration standard 0.080 AC
- Results: 0.086AC & 0.087AC
- Instrument operating properly and accurately. YES NO

COMMENTS: *changed date 07/10/04 → 07/14/04

** moved clock ahead 1020 → 1023

SIGNATURE

Jennifer L. Combs

o. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR
SAMPLE CAPTURE.

P. 5 of 6
flc

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE *b1*

Agency Tucson Police Department Date 7/14/04 Time 1021
Intoxilyzer Serial# 66-005003 Location DUI HQ
QA Specialist Jennifer Combs 461160
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made.
Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Calibration check completed.
Test results 0.086 AC & 0.087 AC
- (✓) 10. Air blank completed.
- (✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Jennifer L. Combs

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

186
P

EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R3-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES ^{A4}
INTOXILYZER MODEL 5000*

AGENCY Tucson Police QA SPECIALIST Irene Doran 46331
(Print Name)
INTOXILYZER SERIAL #46-005003 LOCATION DUE HQ
DATE August 11 2004 TIME 1008

✓ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- ✓ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .394.
- ✓ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- ✓ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "?".
- ✓* Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- ✓ Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- ✓ Alcohol free subject Test result 0.000 AC
- ✓ Error recognition logic system functioning
 Invalid test printed
- ✓ Proper sample recognition system
 Invalid sample printed No
 Deficient sample printed No
- N/A Completeness of sample purge with collection tube**
- ✓ Calibration standard 0.080 AC
- Results: 0.087AC + 0.088AC
- Instrument operating properly and accurately. YES ✓ NO

COMMENTS: #1018 → 1021 moved clock ahead

SIGNATURE

Irene Doran

O. N/A* alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR SAMPLE CAPTURE.

Tucson Police Department Crime Laboratory

396
64

CMI INC
INTOXICLYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-005003
08/11/84 10:43

TEST	BRAC VALUE	TIME
AIR BLANK	.000	10:42
CAL, CHECK	.087	10:42
AIR BLANK	.000	10:42

SUBJECT NAME	DUT THG
NOTE THIS COPYRIGHT	46331
OPERATOR	
ANALYST USE ORGANIC ANALYSIS	
INSTRUMENT PRINTER CARD	✓
INSTRUMENT PRINTER CARD	✓
© 1986 by CMI INC	

CMI

396
64

576
JL

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE AM

Agency Tucson Police Department Date 8/11/04 Time 1008
Intoxilyzer Serial# 66-005003 Location DUI HQ
QA Specialist Irene Doran (46331)
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made.
Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Calibration check completed.
Test results 0.087 + 0.088 AC
- (✓) 10. Air blank completed.
- (✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Irene Doran

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD